

Elite First Aid STOMP Medical Kit



A professional-grade mobile medical station engineered for mass casualty incidents, remote expeditions, and advanced field trauma care — containing over **500 medical items** in a single deployable pack.

🛡️ PROFESSIONAL GRADE

✓ HSA/FSA ELIGIBLE

📁 500+ ITEMS

Mission Briefing & Specifications

The STOMP (Seal Team Operational Medical Pack) transforms a single operator into a capable field medic, providing the tools necessary to manage everything from minor injuries to severe, life-threatening trauma where definitive care is delayed. The kit contains over 500 medical items across 8 categories.

Pack Dimensions & Weight

19" × 14" × 10"

Height × Width ×
Depth

~20 lbs.

Fully Stocked Weight

Construction Details

- **Material:** Heavy-duty rugged nylon construction
- **Zippers:** YKK self-repairing zippers
- **MOLLE Webbing:** 5 external rows for modular expansion
- **Carry System:** Padded adjustable backpack straps with sternum & waist support
- **Handles:** Dual top carry handles for rapid grab-and-go


Bleeding Control & Trauma Loadout

The hemorrhage control suite is the cornerstone of this kit — engineered to stop life-threatening bleeding fast with professional-grade hemostatic agents and pressure dressings.



Hemostatic Agents

- 3 × QuikClot 4"x4" Dressing
- 4 × BleedStop Bandages
- 2 × WoundSeal Powder Packs

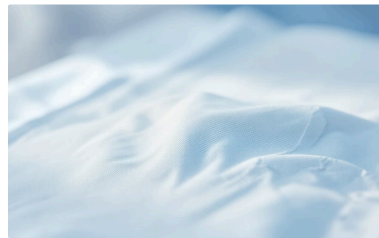

Pressure Dressings

- 2 × 4" Trauma Dressing
- 2 × 6" Trauma Dressing
- 2 × 12"x30" Multi-Trauma Dressing



Hemorrhage Control

- 10 × 5"x9" Trauma Pads

Full Component Loadout — 500+ Items

Every category of medical need is addressed, from airway management and surgical intervention to burns, fractures, and general wound care.

Bandages & Gauze

- 100 × 1"x3" Bandages
- 20 × 2"x4" Bandages
- 20 × Knuckle Bandages
- 20 × Butterfly Closure Bandages
- 5 × Wound Closure Strips
- 20 × 2"x2" Gauze
- 10 × 4"x4" Gauze
- 2 × 4.5"x4yd Krinkle Gauze
- 3 × Triangular Bandages
- 5 × Eye Pads

Airway & CPR

- 2 × CPR Face Shields
- 2 × 80mm Oral Airways
- 2 × 110mm Oral Airways

Diagnostics & Vitals

- 1 × Blood Pressure Unit
- 1 × Stethoscope
- 1 × Penlight
- 5 × Tongue Depressors

Instruments & Surgical Tools

- 2 × Trauma Shears
- 1 × Tweezers
- 2 × Straight Hemostats
- 2 × Curved Hemostats
- 1 × Scalpel Handle
- 6 × Scalpel Blades
- 1 × Needle Probe

Orthopedic, Splinting & Wraps

- 1 × Cervical Collar
- 2 × Splints
- 6 × 2"x5yd Elastic Bandages
- 3 × 6"x5yd Elastic Bandages
- 3 × 1"x10yd Tape
- 1 × Elastic Strap
- 3 × Cold Packs

OTC, Topicals & Burn Care

- 20 × Ibuprofen 2/packs
- 20 × Triple Antibiotic Packets
- 1 × Burn Spray
- 5 × Burn Gel
- 2 × Anti-Itch Packets
- 12 × Sting Wipes
- 20 × Ammonia Wipes
- 1 × Glucose Pack
- 2 × Lip Balm

Antiseptics, Cleaning & Eye Care

- 100 × Iodine Wipes
- 28 × BZK Antiseptic Wipes
- 20 × Alcohol Wipes
- 1 × 2oz Hand Sanitizer
- 2 × Eyewash

Personal Protection, Survival & Misc.

- 5 pair × Gloves
- 2 × Emergency Blankets
- 6 × Safety Pins
- 1 × First Aid Guide

500+

Total Items

Complete medical loadout
in one pack

8

Supply Categories

Bandages, Trauma,
Antiseptics, OTC/Burn,
Instruments, Diagnostics,
Airway, Orthopedic

3

QuikClot Dressings

Hemostatic dressings
included in the kit

100

Iodine Wipes

Antiseptic coverage for
mass casualty events

Tactical Deployment & Features

Every design decision in the STOMP pack prioritizes speed, organization, and operational effectiveness under pressure.

Rapid Access Tri-Fold Design

The pack unzips completely to lay flat, exposing the **entire loadout simultaneously**. This eliminates digging for supplies, allowing rapid visual identification and access to every tool via mesh pockets, elastic loops, and removable pouches.

Advanced Trauma Capabilities

Multiple hemostatic agents, Israeli Trauma Dressings, and a tourniquet provide immediate hemorrhage control. The surgical set, sutures, SAM splints, and cervical collar allow an operator to **stabilize fractures, close wounds, and manage complex injuries** in the field.

Modular Organization & Expansion

Two main interior sections feature mesh pockets, vertical slide pockets, instrument panels, and removable zippered pouches. Externally, **5 rows of MOLLE webbing** allow attachment of additional gear, while 6 adjustable compression straps secure contents during transport.

Who Carries the STOMP?

The Elite First Aid STOMP is purpose-built for operators and responders who cannot afford to be under-equipped when lives are on the line.



Tactical Medics & EMS

Professional first responders who need a comprehensive mobile medical station deployable in seconds at any incident scene.



Remote Expedition Leaders

Guides and expedition leaders operating in wilderness or off-grid environments where medical help may be hours or days away.



Disaster Relief & SAR Teams

Search and rescue operators working in chaotic, resource-scarce environments where the kit must serve as a complete field clinic.



CERT & Preppers

Community Emergency Response Teams and serious preppers who understand that comprehensive preparation is the foundation of resilience.

Field Scenarios & Recommended Environments

The STOMP is designed to perform across the full spectrum of emergency environments — from urban mass casualty events to remote wilderness operations.

Mass Casualty Incidents

Provides the breadth of supplies needed to treat **multiple patients** simultaneously. The 500+ item loadout ensures you won't run dry in a multi-victim scenario.

Remote & Austere Environments

Serves as a **mobile clinic** for expeditions, hunting parties, or off-grid living where medical help is hours or days away. Every critical intervention is covered.

Natural Disaster Response

An essential component for any go-bag or vehicle kit during **hurricanes, earthquakes**, or other emergencies. Grab it and go when seconds count.

Tactical Operations

Provides comprehensive medical support for **law enforcement or military units** operating in high-risk environments where trauma care must be immediate.

Maintenance & Re-Stowage Protocol

A medical kit is only as reliable as its maintenance. Follow this protocol to ensure the STOMP is always mission-ready.

Quarterly Inspection Checklist

01

Inventory Check

Lay out all contents and verify against the component list. Order replacements for any used or missing items immediately.

02

Expiration Check

Systematically check expiration dates on all dated materials — hemostatic agents, medications, and sterile dressings. Rotate or replace items nearing expiration.

03

Damage Assessment

Inspect the backpack for tears, broken buckles, or failing zippers. Check all tools and instruments for rust or damage.

Re-Stowage Procedure

01

Prioritize Access

Repack with the most critical life-saving items — tourniquets, pressure dressings, hemostatics — in the most accessible locations.

02

Logical Grouping

Keep like items together in designated pouches: all airway supplies in one area, all diagnostic tools in another.

03

Secure Contents


Utilize internal cinch straps to secure larger items. Ensure all pouches are properly zipped and secured before closing.

04

Final Cinch

After zipping the pack closed, tighten all **6 external compression straps** to secure the load and minimize profile for transport.

Maintenance Timeline at a Glance

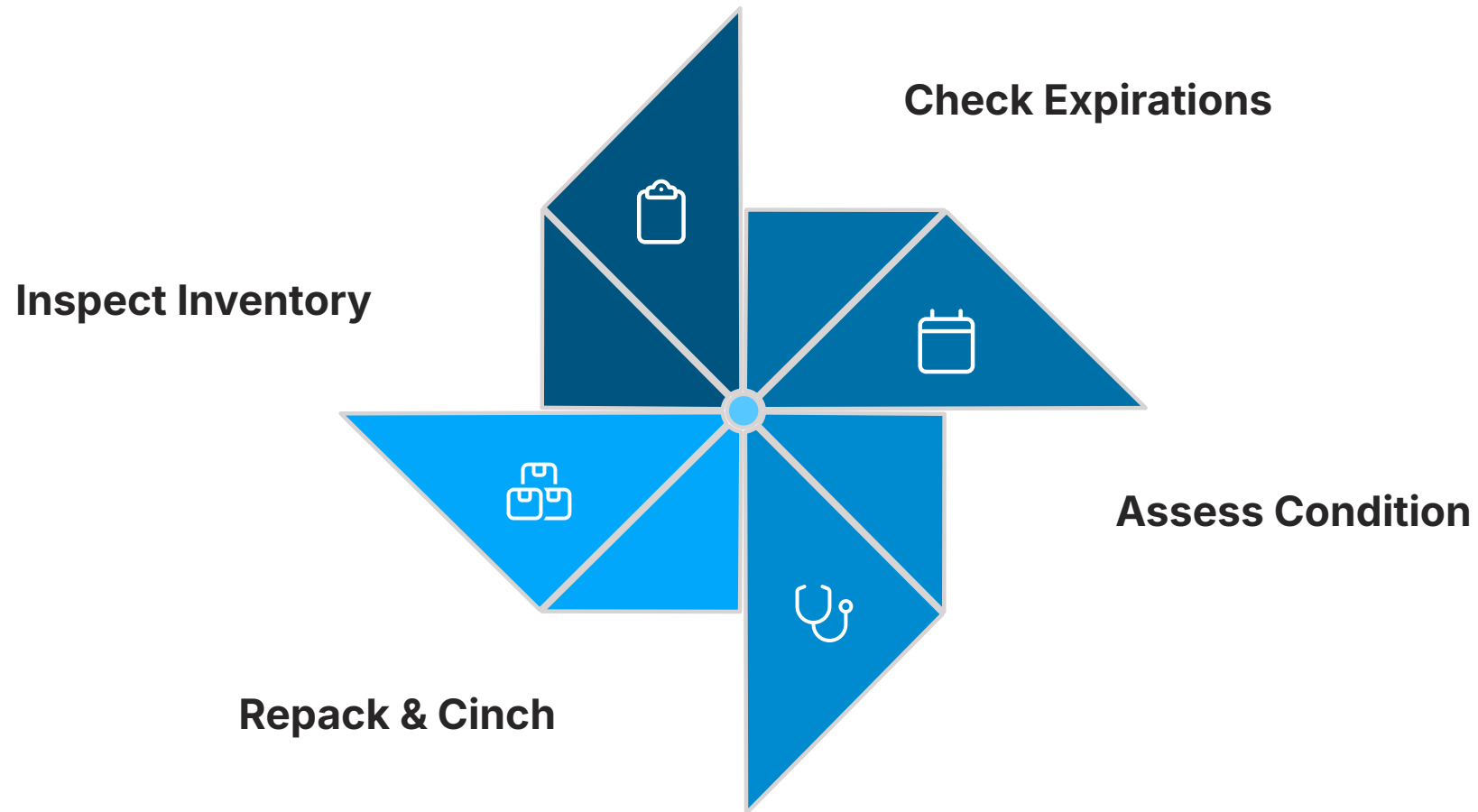

Consistent quarterly maintenance ensures the STOMP is always mission-ready. A neglected kit is a liability — a maintained kit is a lifeline.

- Pro Tip:** Mark your quarterly inspection dates on a calendar and store a laminated copy of the component checklist inside the front pocket of the pack. This ensures any operator — not just the primary owner — can conduct a proper inventory check in the field.

The STOMP — Mission Ready, Every Time

The Elite First Aid STOMP Medical Kit is more than a first aid bag. It is a **deployable medical station** that empowers a single trained operator to deliver advanced trauma care when and where it matters most.

500+ Items

Comprehensive loadout across 8 medical categories

Tri-Fold Access

Full loadout visible and accessible in seconds

MOLLE Ready

5 rows of external webbing for modular expansion

HSA/FSA Eligible

Qualified medical expense for tax-advantaged purchase

Engineered for the operators who understand that in a crisis, the quality of your equipment is the difference between life and death. The STOMP is built to that standard — without compromise.

